

Synonym

ACVR2B,ACTRIIB,MGC116908

Source

Human Activin RIIB, His Tag(ACB-H5226) is expressed from human 293 cells (HEK293). It contains AA Ser 19 - Thr 137 (Accession # [Q13705-1](#)).

Predicted N-terminus: Ser 19

Molecular Characterization

ACVR2B(Ser 19 - Thr 137)  
Q13705-1

Poly-his

This protein carries a polyhistidine tag at the C-terminus.

The protein has a calculated MW of 14.5 kDa. The protein migrates as 20-35 kDa when calibrated against [Star Ribbon Pre-stained Protein Marker](#) under reducing (R) condition (SDS-PAGE) due to glycosylation.

Endotoxin

Less than 1.0 EU per µg by the LAL method / rFC method.

Purity

>90% as determined by SDS-PAGE.

>90% as determined by SEC-MALS.

Formulation

Lyophilized from 0.22 µm filtered solution in PBS, pH7.4 with trehalose as protectant.

Contact us for customized product form or formulation.

Reconstitution

Please see Certificate of Analysis for specific instructions.

*For best performance, we strongly recommend you to follow the reconstitution protocol provided in the CoA.*

Storage


For long term storage, the product should be stored at lyophilized state at -20°C or lower.

*Please avoid repeated freeze-thaw cycles.*

This product is stable after storage at:

- 20°C to -70°C for 12 months in lyophilized state;
- 70°C for 3 months under sterile conditions after reconstitution.


SDS-PAGE

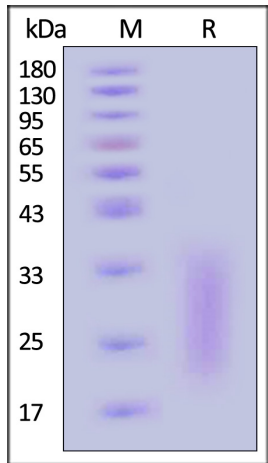

Human Activin RIIB, His Tag on SDS-PAGE under reducing (R) condition. The gel was stained with Coomassie Blue. The purity of the protein is greater than 90% (With [Star Ribbon Pre-stained Protein Marker](#)).

Bioactivity-SPR

SEC-MALS

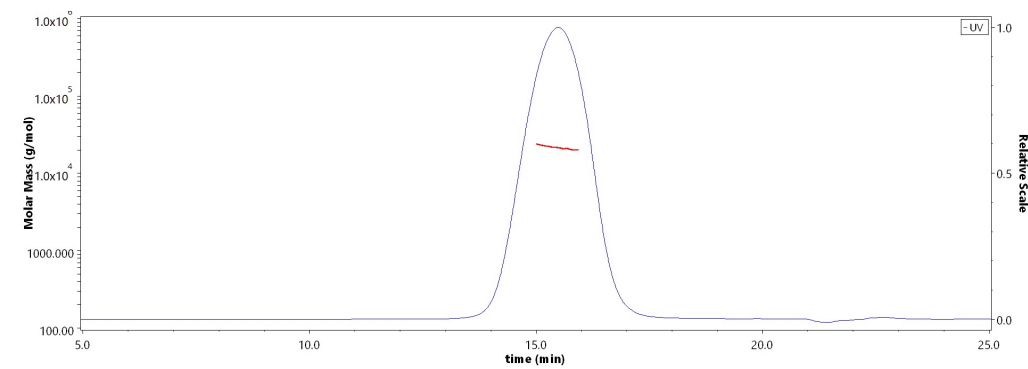


The purity of Human Activin RIIB, His Tag (Cat. No. ACB-H5226) is more than 90% and the molecular weight of this protein is around 18-25 kDa verified by SEC-MALS.

[Report](#)

Discounts, Gifts,  
and more!



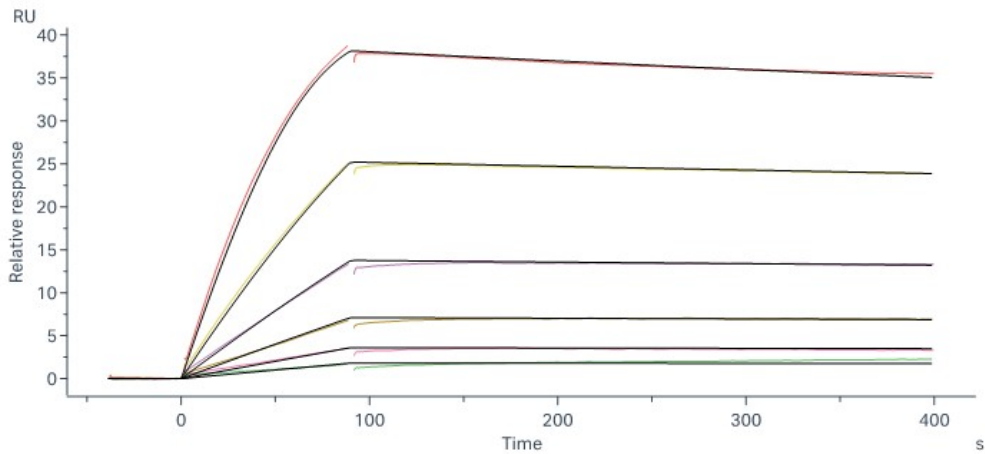

Human Activin RIIB, His Tag (Cat. No. ACB-H5226) immobilized on CM5 Chip can bind Human Activin A, premium grade (Cat. No. ACA-H421b) with an affinity constant of 0.164 nM as determined in a SPR assay (Biacore 8K) (QC tested).

Background

Activin receptor type-2B (ACVR2B) is also known as ActR-IIB and MGC116908, ACVR2B is an activin type 2 receptor. Activins are dimeric growth and differentiation factors which belong to the transforming growth factor-beta (TGF-beta) superfamily of structurally related signaling proteins. Activins signal through a heteromeric complex of receptor serine kinases which include at least two type I (I and IB) and two type II (II and IIB) receptors. These receptors are all transmembrane proteins, composed of a ligand-binding extracellular domain with cysteine-rich region, a transmembrane domain, and a cytoplasmic domain with predicted serine/threonine specificity. Type I receptors are essential for signaling; and type II receptors are required for binding ligands and for expression of type I receptors. Type I and II receptors form a stable complex after ligand binding, resulting in phosphorylation of type I receptors by type II receptors. Type II receptors are considered to be constitutively active kinases. This gene encodes activin A type IIB receptor, which displays a 3- to 4-fold higher affinity for the ligand than activin A type II receptor. Defects in ACVR2B are the cause of visceral heterotaxy autosomal type 4 (HTX4).

